


- Framgaffelfästet, justerbart i höjd och vinkel.
- Byte av länkhjulets läge i framgaffeln.
- Byte av framgaffel
- Byte av länkhjul
- Byte av framgaffelfäste

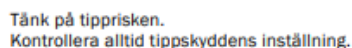

			Länkhjulsposition (cm)											
			Gaffel extra kort			Gaffel Kort			Gaffel Medium			Gaffel Lång		
			1	2	3	1	2	3	1	2	3	1	2	3
Lågt framgaffelfäste	H	4" *	32-36	33.5-37.5	-	-	-	-	-	-	-	-	-	-
	G	4" *	- 34-38	35.5-39.5	-	-	-	-	-	-	-	-	-	-
	G	5"	- -	-	36-40	37.5-41.5	-	-	-	-	-	-	-	-
	A-F	5"	- -	-	36-40	37.5-41.5	39-43	-	-	-	-	-	-	-
		6½"	- -	-	-	-	41-45	-	-	-	-	-	-	-
Standard/ Breddat framgaffelfäste	A-F	5"	- -	-	37.5-41.5	39-43	40.5-44.5	43-46	45-48	46-49	-	-	-	-
		6"	- -	-	-	40-44	41.5-45.5	45-47	46-49	47-50	-	-	-	-
		6½"	- -	-	-	41-45	42.5-46.5	46-48	47-49	48-51	50.5-52.5	51.5-53.5	52.5-55.5	-
	A-E	7"	- -	-	-	-	-	-	47.5-49.5	48.5-51.5	-	-	-	53-56
	A-C	8" **	- -	-	-	-	-	-	48-51	49-52	-	52.5-54.5	53.5-56.5	-

För länkhiul 8x2", använd gaffel Transit (ger sitshöjd 47-51 cm).

Bakre sitshöjd kan justeras med hjälp av:

- Drivhulets (camberbrickans) position i gaveln.
- Camberbrickans vridning i gaveln.
- Byte av drivhjul




Tänk på tipprisken.
Kontrollera alltid tippskyddens inställning.

Bakre sitshöjdstabell			
	Drivhjul	Sitshöjd cm	Avstånd cm sits-drivhjul
Position A	20"	46	4.5
	22"	48.5	7.0
	24"	51	9.5
	26"	53,5	12
Position B	20"	44	6.5
	22"	46.5	9.0
	24"	49	11.5
	26"	51,5	14
Position C	20"	42	8.5
	22"	44,5	11.0
	24"	47	13.5
	26"	49,5	16
Position D	20"	40	10.5
	22"	42.5	13
	24"	45	15.5
	26"	47,5	18
Position E	20"	38	12.5
	22"	40.5	15
	24"	43	17.5
	26"	45,5	20
Position F	20"	36	14.5
	22"	38.5	17
	24"	41	19.5
	26"	43,5	22
Position G	20"	34	16.5
	22"	36.5	19
	24"	39	21.5
	26"	41,5	24
Position H	20"	32	18.5
	22"	34.5	21
	24"	37	23.5
	26"	39,5	26



Sitshöjd

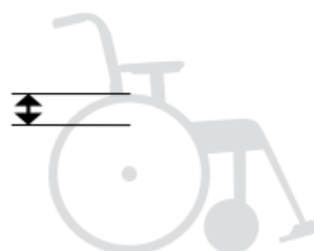

Avstånd sits-drivhjul